



ECE Group II - III
15 - 36 kg
33.1 - 79.4 lb

COMMANDER

The Specialist – Highly Adaptive Positioning Car Seat for children with special needs



THOMASHILFEN COMMANDER

” The Thomashilfen Commander masters every challenge.
As a special seat for special requirements and positioning, it is ideally suited for highly individual supply solutions. Children are supported wherever they need support and protection. Especially for children with little upper body control, the fixed lateral trunk supports offer this extra support – the supports have proven themselves over many years of use in rehab strollers. A huge selection of accessories and easy handling make the adjustment perfect. “



Particularly suitable for:

- Cerebral Palsy (CP)
- Skull-brain trauma up to grade II
- Spina Bifida
- Rett Syndrome
- Autism
- Spinal Muscular Atrophy (SMA)

Always recommended for


- Children with little upper body control



Thomashilfen has been providing children with special needs with rehab car seats for over 30 years. Today we also use our know-how to develop our own rehab car seats. All our many years of experience with the needs of children and parents flow into them.

A handwritten signature in black ink that reads "Gunnar Thomas".

Gunnar Thomas, owner

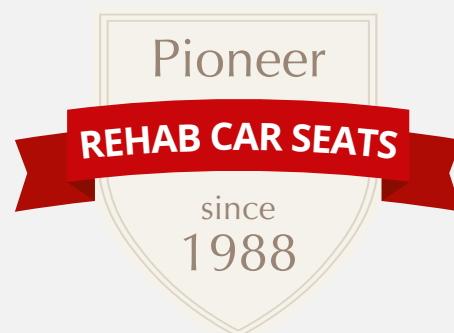

REHAB-EQUIPMENT



UNIQUE COMBINATION



The unique, modular combination of **SWIVEL BASE** and **TILTING BASE** can be supplemented with a variety of accessories as required.



For easy handling, the release lever (draw latch) on the swivel base is easy to reach. It can be mounted either on the left or right front of the swivel base.

UPHOLSTERY

cuddly upholstery for pleasant sitting comfort without annoying sweating, the breathable cover can be completely removed and washed at 40° C



CUP HOLDER

removable large holder (right or left) to keep drinks, snacks and toys ready to hand



SEAT TILT-IN-SPACE

with the integrated seat tilt-in-space +3° a resting position is quickly adjusted

HEAD PADS

additional padding inserts can be added on the right and left side of the head pads, even better support for the child's head



HEAD SUPPORT

9-way height adjustable – easy and quick to adjust, grows step by step with your child



SHOULDER BELT

"grows" with the child, without the need to rethread the belts



BACK PLATE

sturdy, large back plate with hook-and-loop-fasteners for attaching the upper body supports – for optimum individual adjustment

Two versions **LATERAL TRUNK SUPPORTS** to choose from:



FOAM CORE: soft lateral trunk pads to attach directly onto the cover with hook-and-loop-fasteners



METAL CORE:

stable lateral trunk supports for attaching to the back plate with hook-and-loop-fasteners

5-POINT-BELT

integrated 5-point positioning belt with central adjuster for quick adjustment of the belt length, practical chest clip prevents twisting of the belts



ACCESSORIES – COMBINATION



ADAPTER PLATE –
ALWAYS required for the
following accessories

The unique modular combination
of swivel base and tilting base can
be combined with a variety of
accessories.

**An adapter plate is ALWAYS
required to use the accessories.**



SWIVEL BASE facilitates the
transfer, easily accessible
draw latch can be used left
or right



TILTING BASE allows seat
tilt in 3 steps (5.5°, 11°,
16°) for relaxed resting
position



TILTING BASE and **SWIVEL
BASE** combined for best
possible handling and
comfort



FOOTREST WITH PAD
provides stable sitting posi-
tion, children can firmly
support their feet and do
not slip out of the seat



STAND, SEATFIX ADAPTER
and **STABILIZING BOW**
guarantee stable seat
support in the car



SEAT DEPTH EXTENSION
allows the seat to grow
with the child and ensures
correct positioning of the
thighs



ABDUCTION BLOCK
(adjustable) ensures stable
sitting position and pre-
vents legs from crossing
over

THOMASHILFEN Commander



All details online:
www.thomashilfen.com/commander



Accessories



ABDUCTION BLOCK and **HIP SUPPORTS** support a stable sitting position



LATERAL TRUNK SUPPORTS (METAL CORE) can be positioned precisely on the back plate with hook-and-loop-fasteners for individual adjustment



CROTCH PAD (LARGE) protects the sensitive belly of the child, **SEAT WEDGE (BELOW)** allows +10° seat tilt



SEAT WEDGE (INSIDE) supports an upright pelvic position, **LATERAL TRUNK SUPPORTS (METAL CORE)** for optimum support of the upper body,

Accessories



TABLE promotes the straightening up and stable posture of the child, provides the child with all-round support and safety



ANTI-ESCAPE VEST prevents the child from wriggling out of the belt



CAR SEAT PROTECTION MAT WITH BACK protects the seat upholstery and offers 3 practical compartments for toys at the front

A complete overview of accessories can be found on the order form on our website:
www.thomashilfen.com

Dimensions & weights	THOMASHILFEN Commander
Seat depth	34 cm / 13.4"
Seat width front (with supports)	33,5 cm (22 cm) / 13.2" (8.7")
Seat width back (with supports)	28,5 cm (17 cm) / 11.2" (6.7")
Back height (seat)	55 - 75 cm / 21.7 - 29.5"
Shoulder height (belt guide)	25 - 44 cm / 9.8 - 17.3"
Seat tilt-in-space (basic model)	3°
(with tilting plate)	5,5° / 11° / 16°
Lower leg length	23 - 40 cm / 9.1 - 15.7"
(Footrest directly attached to adapter plate)	
Lower leg length	28 - 45 cm / 11 - 17.7"
(Footrest attached to swivel base + tilting base)	
Overall dimensions (w x h x d)	55 x 65 x 45 cm / 21.7 x 25.6 x 17.7"
Weight of seat	approx. 8,5 kg / 18.7 lb
max. user weight	15 - 36 kg / 33.1 - 79.4 lb



5 upholstery designs available:



Basic configuration:

Seat and back with integrated 9-fold height-adjustable headrest; back plate (with hook-and-loop fasteners); seat cushion; back cushion; head cushion; additional head cushion inserts; crotch pad, small; 5-point positioning belt with central adjuster; chest clip; shoulder protector; cup holder (can be used right or left); removable covers, machine washable up to 40° C



Thomas Hilfen für Körperbehinderte GmbH & Co. Medico KG · Walkmühlenstraße 1 · D - 27432 Bremervörde
 Phone: +49 (0)4761 886-63 or -68 · Fax: +49 (0)4761 886-19 · www.thomashilfen.com